



Composite placement and modelling instrument


Compothixo

Compothixo – Improved Quality Restorations

- Better wettability
- Superior adaptation of composite to cavity walls
- Reduction of air bubbles
- Precise application
- Layer thickness control
- Improved sculptability
- Reduced stickiness

Compothixo represents a unique generation composite placement and modelling instrument suitable for all classes of restorations.

New **Compothixo** technology enhances the thixotropic properties of composites by only changing its viscosity, without altering the chemical and mechanical characteristics of the material.

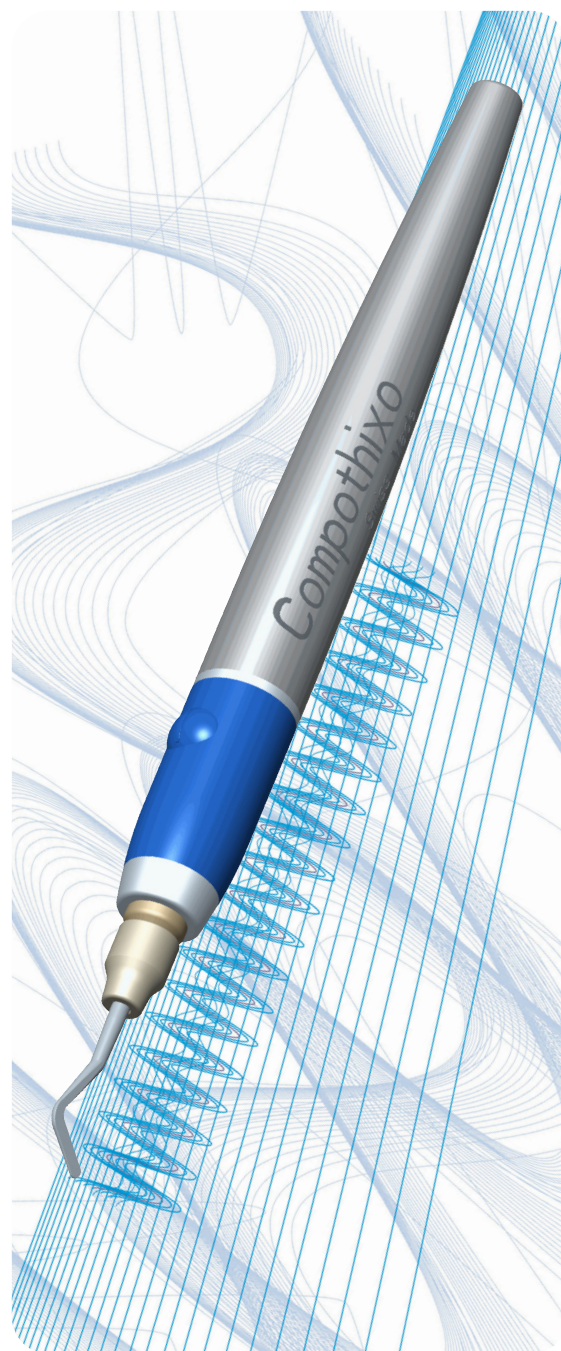

Your practice is our inspiration.™



Composite placement and modelling instrument

Compothixo

Indications:

- Modelling of composite
- Occlusal modelling, fissures and removal of excess
- Layer application technique
- Bulk technique in small cavities
- Direct veneering

Order Information

5400 Compothixo Intro Kit

Contains: Instrument handle + 4 Tips (1 Pointed Tip, 1 Spatula Tip, 1 Plunger Tip, 1 Semi-Sphere Tip)

5401 Compothixo Instrument Handle

Contains: 1 Instrument handle

5402 Compothixo Tips Assorted Kit

Contains: 1 Pointed Tip, 1 Spatula Tip, 1 Plunger Tip, 1 Semi-Sphere Tip

Refills Tips:

Contains: 2 pcs

5403 Compothixo Pointed Tip

5404 Compothixo Spatula Tip

5405 Compothixo Plunger Tip


5406 Compothixo Semi-Sphere Tip

Benefits :

- Superior cavity and matrix adaptation
- Reduction of air bubbles
- Easy thin layer stratification, especially in anterior restorations
- Reduced pressure required for packing of composite
- Reduced stickiness of the composite versus traditional instruments
- Easy composite placement on the tooth wall; no pull-back effect
- Easier application of bulk technique in small cavities

Technical features:

- Ergonomic design
- Light weight
- Four common used interchangeable tips
- Easy tip fitting
- Battery driven
- 140 Hz vibration frequency
- 8 hours battery autonomy

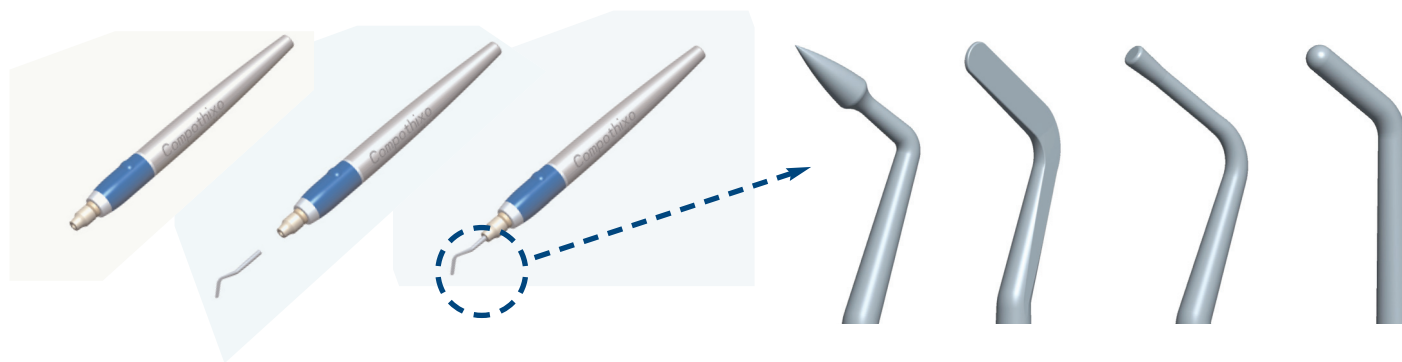

KerrHawe SA

Via Strecce 4, P.O. Box 268, 6934 Bioggio, Switzerland

Freephone: 00800 41 05 05 05 Fax: ++41 91 610 05 14 www.kerrdental.eu

Kerr
Sybron Dental Specialties